

# Bariatric Range

Rental Range and Options

permobil



# About Permobil

At Permobil we put the people who use our products first. Our purpose is to innovate for individuals; to create and provide advanced assistive solutions that make the lives of people living with disabilities more enriching. Our dedicated teams work tirelessly to make this happen, no matter if they are designing a new wheelchair, testing seating cushions, or supporting a therapist.

Our founders were passionate about providing a better life for people with mobility impairments and that same passion still flows through everyone at Permobil. It drives us to push the boundaries of technology and our own ingenuity in pursuit of innovations that improve the quality of life of our users. Just as our users' needs constantly change, so we work continuously to tailor Permobil products to meet their day-to-day reality.

Thanks to this, our users can trust that we provide market leading solutions, driven forward by working with insights based on clinical data.

At Permobil, we are leading change. We see ability. We are innovating for individuals.

This catalogue outlines the range we currently offer to support bariatric users. We have included products that cater for users who require a SWL over 180kg. The exception is paediatric equipment and specific information is provided with those items.

## Our core values guide us in our response to our customers needs:



**Tuatahi**  
User First



**Hautūtanga**  
Leading










**Whakawhirinaki**  
Trusted



**Ngākau Whiwhita**  
Passionate

# Contents

About Permobil	2
 Bed Range	4
 Mattress Range	8
 Wheelchair Range	12
 Cushion Range	14
 Chair Range	16
 Bathroom Range	18
 Manual Handling Range	21
Permobil Clinical Resources	23
FAQs and Contacts	24



Represents PHARMAC contracted items. Visit their website by scanning:



## Contact Us

### Auckland Rental

Rental.NZ@permobil.com

0800 115 222

181A McLeod Road, Te Atatu South

Auckland 0610, New Zealand

### Wellington Rental

RentalNZ.Nth@permobil.com

0800 363 123

15 Raroa Road, Hutt Central

Lower Hutt 5010, New Zealand

### Christchurch Rental

RentalNZ.Sth@permobil.com

03 354 9239

3/114 Sawyers Arms Road, Northcote

Christchurch 8052, New Zealand

### Permobil Contract Manager

Susie Ballantyne

Commercial Services Director

Susie.Ballantyne@permobil.com

021 246 1885

# Bed Range

We have a range of beds available in single or king single sizes.

## Bariatric Floorline Bed - Big Ted

Code: BED BARI RENTAL **P**

The Big Ted Bed is a three function bed with a generous mattress platform, consisting of four panels.

- Fitted with twin actuators on back rest and knee bend functions, with pair of strong actuators at both ends of bed to raise/lower the mattress platform.
- Pump bracket is fitted to foot end tower to accommodate an alternating mattress.
- The bed can be folded for ease of transport.
- Available in King Single.



Specifications	
SWL	350kg
Width	107cm
Length	203cm
Height	13.5cm to 73.5cm
Profiling Functions	3

## Botero Bariatric Bed

Code: BOTERO BARI BED **P**

The specific Bariatric design offers a width and length adjustable mattress platform to manage changing needs in both ward and respite environments.


- Low height of 38cm.
- Attendant control panel with intuitive pre-programmed single touch functions for speed.
- Central locking castors with integrated directional lock.
- Available in King Single.
- Can extend to Super King Single as width is adjustable.



Specifications	
SWL	400kg
Width	102cm - 122cm
Length	226cm - 244cm
Height	38cm - 85cm
Trendelenburg	+17 degrees/ -17 degrees
Profiling Functions	4

# Bed Range

## Floorline Bed - Rose II

Code: BED S FLOOR  
BED KS 

The Rose II allows patients to be lowered to floor level and raised to a suitable height.

- Electrically operated backrest, knee break, raise/lower and tilt functions.
- Footend tower has built in keypad to operate bed and separate hand control comes standard.
- Bed can be folded for storage, transport and installation and fits through standard doorways.
- Available in Single or King Single.



Specifications	
SWL	220kg
Single	W 90cm x L 193cm
King Single	W 107cm x L 203cm
Height	12.5cm - 72cm
Profiling Functions	4

## AC3 Bed/AC4 Bed

Code: BED AC3  BED AC4 KS

Both beds are designed with the care of residents and caregivers in mind, focusing on ease of use, functionality, and reliability.

- Simple and intuitive.
- Added comfort of a larger sleeping surface.
- Quick attaching and easily removable head and foot boards.
- AC3 - Available in Single.
- AC4 - Available in King Single.



Specifications	
SWL	220kg
Single	W 92cm x L 210cm
King Single	W 107cm x L 210cm
Height	30cm - 73cm

# Bed Range

## Walmsley Bed

Code: BED S WALMSLEY

The Walmsley range of healthcare beds are a commercial grade, electrically operated hi-low bed, which have been designed for maximum patient /carer safety and reliability.

- The electrical components are contained under the mattress platform for safety and protection while all critical moving parts are fitted with bronze bearings for a quiet and smooth operation.
- Wheels have both brake and directional lock, once the locking tab is engaged.
- An easy to operate coloured hand control is standard on all Walmsley beds.
- Easy to remove head and foot boards are standard.
- Available in Single.



Specifications	
SWL	250kg
Width	92cm
Length	196cm
Height	24cm - 74cm

## Bed Accessories

Bed accessories are available on some of our rental beds. If you need any of the below, please reach out and ask.

Code: - BED STICK RENTAL    Bed stick **(P)**

Code: - BED RAILS RENTAL    Bed side rails - single or pair **(P)**

Code: - SELF HELP RENTAL    Self-help pole **(P)**

Please note we can also supply IV poles and Oxygen bottle holders for some beds.

# Bed Range

## ook Snow ALL Bariatric Hospital Bed

Code: BED OOK

Featuring the ability to expand in width, institutions can meet the needs of bariatric, medical-surgical, and palliative care settings. With seamless smooth surfaces and liftable deck covers, the ook snow ALL bariatric hospital bed has been thoughtfully designed to maximise infection control. One-touch buttons provide ease of use. The exceptional low height, bed exit alarms, integrated nurse call possibilities and scales support help reduce the risk of falls and improve patient monitoring.

- Sleeping surface dimensions: W 89cm to 119.4cm, L 203cm to 213cm
- Trendelenburg: -16° to +16°
- Backrest angle: 0 to 60°
- Smartscreen: Optimise patient safety through advanced bed monitoring capabilities
- Low height: Exceptional low height reduces the risk of falls and facilitates patient mobility
- Easy to clean: Open architecture, seamless smooth surfaces, and liftable deck covers facilitate cleaning procedures



Specifications	
SWL	481kg
Width	105.4cm to 134cm
Length	228.6cm to 238.7cm
Height	28cm - Low Height 76cm - High Height
Profiling Functions	4

## ook Snow ALL Bariatric Hospital Bed with Pilot Drive

Code: BED OOK PILOT

As above but with a powered drive system designed with mobility, versatility, and patient safety all-in-one.

## Bed Accessories

Code: OOK HELPER





Code: OOK PADS SET



# Mattress Range

Our mattresses are suitable for a range of support requirements, bed sizes and care environments.

## Active Air BARI 10 Mattress Replacement

Code: ASPIRE BARI S, ASPIRE BARI SK   
ASPIRE BARI SKS 

- Adaptable Pump - can be switched between modes for users up to 160kg, and bariatric users up to 450kg.
- Integrated Turning Therapy - allows patient rotation up to 40° for improved pressure care, assists staff with patient handling requirements.
- Rapid Pump Inflation - maximises therapy response times, minimises staff resources by fully inflating in only two minutes.
- Hospital Grade Cover - thicker weight polyurethane fabric with increased bleach resistance to handle modern hospital infection control standards.
- Available in Single, King Single or Super King Single.



Specifications	
SWL	450kg
Single	90cm x 200cm
King Single	107cm x 200cm
Super King Single	123cm x 200cm
Depth	27cm
Cell Cycle	1 in 2, 1 in 4
Cycle Time	10 / 15 / 20 mins

## Forté Airmonte

Code: MA AIRMONTE S

The Airmonte is easy to use and comfortable, making it an ideal solution for any care environment aiding in the prevention and treatment of pressure injuries.

- Very quiet pump operation.
- Auto-Adjust: Automatically adjusts pressure to suit patient's weight.
- Auto-Firm: Switches to static mode and max-inflates to allow for easy repositioning and transfers.
- Light overall weight (9kg) allows ease of set-up and transportation.
- Available in Single.



Specifications	
SWL	30kg - 230kg
Width	87cm
Length	198cm
Depth	20cm
Cell Cycle	1 in 2
Cycle Time	10 / 15 / 20 / 25 mins

# Mattress Range

## QUATTRO® Acute Mattress System

Code: ACUTE **(P)** ACUTE KS **(P)**

- Features include Active QUATTRO THERAPY and Continuous Low Pressure Therapy modes, together with AUTOsupport fully automatic pressure adjustment and the unprecedented low interface pressures with DEEP CELL THERAPY.
- Available in Single or King Single.



Specifications	
SWL	250kg
Size	Dimensions (W x L x D cm)
Single	90cm x 198cm x 24cm
King Single	110cm x 198cm x 24cm
Number of Cells	27
Alternating Cycle	16 minutes
Active Therapy	1-in-4 Active Cell Cycle

## Forté Alaise

Code: MA ALAISE S, MA ALAISE KS

The Alaise is an advanced Hybrid Static Care Support Surface combining together the pressure redistribution benefits of both foam and air, ideal for long-term care applications where maximum pressure care is required with virtually zero maintenance.

- 4 x Anatomically zoned air cell structure - Scapular, Pelvic, Thigh and Heel Zones.
- Premium firm foundation layer and hinged strengthened sides.
- A non-powered, self-regulating system - minimal or no training required for set up.
- High open-cell elasticity for temperate stability and to maximise functionality and patient mobility.
- Available in Single or King Single.



Specifications	
SWL	30kg - 300kg
Size	Dimensions (W x L x D cm)
Single	88cm x 198cm x 18.5cm
King Single	105cm x 198cm x 18.5cm


# Mattress Range

## Forté Sovereign ROHO®

Code: MA SOV ROHO S, MA SOV ROHO KS

The Sovereign High Care mattress with integrated ROHO® DRY FLOATATION® Technology provides trusted and proven care for high to very high-risk clients requiring a static support surface that balances immersion, breathability, postural support, superior comfort, and outstanding performance.

- High-temperature stability to deliver consistent immersion, performance and comfort.
- ROHO® valve extensions to allow pressure to be adjusted without opening cover resulting in minimising infection risk.
- Firm foundation layer with sloped heel pressure transfer zone and hinged strengthened sides.
- High-Performance, 4 zoned immersion region.
- The number of ROHO® inserts can be adjusted between 1-4 depending on the patient needs.

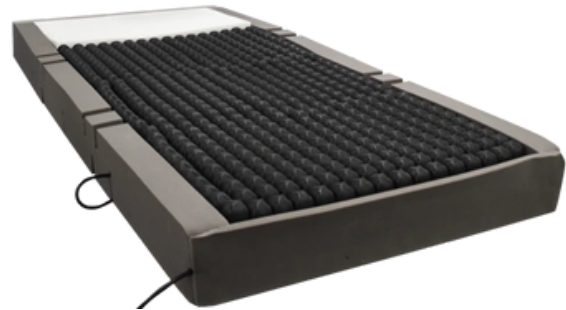


Specifications	
SWL	30kg - 275kg
Size	Dimensions (W x L x D cm)
Single	88cm x 198cm x 15cm
King Single	105cm x 198cm x 15cm

## Modular Therapy Mattress (MTM)

Code: MTM S **P** MTM KS **P** MTM Q **P**

Full mattress replacement using a combination of foam and ROHO® overlay sections.

- Modular system comprised of four overlay sections and foam base. The modular system allows a combination of ROHO® overlay and high-quality foam sections to be used as required for comfort and stability.
- The ROHO® overlay offers excellent pressure redistribution and minimises risk of shearing, while the foam base allows a 12.5cm overall height, similar to a standard mattress.
- The MTM mattress offers a REACTIVE surface and does not require a pump to maintain inflation.
- Available in Single, King Single or Queen.

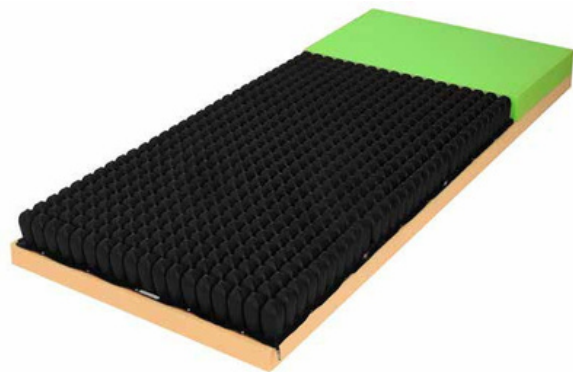

Specifications	
SWL	275kg
Size	Dimensions (W x L x D cm)
Single	87cm x 198cm x 12.5cm
King Single	107cm x 202cm x 12.5cm
Queen	153cm x 202cm x 12.5cm

# Mattress Range

## ROHO® Bariatric Overlay Kit

Code: OVERLAY BARI ROHO 

- No individual weight limit as long as assembled support surface sections are sized to fit the bariatric mattress and the individual.
- ROHO® DRY FLOATATION® overlay cell sections extend width of ROHO® DRY FLOATATION® overlay to fit bariatric bed sizes.
- The four sections of overlay are independently adjusted, providing low pressure environment over the entire contact area of body - available in 3, 4, or 5 cells wide.
- Add on sections help provide the bariatric individual with therapeutic benefits of ROHO® DRY FLOATATION® Technology.



# Wheelchair Range

We have a selection of power and manual chairs to choose from, including Permobil options.

## Permobil M300HD Corpus

Code: PWC M300HD 

- The M300 base is highly maneuverable. The suspension has been strengthened and base extended for a more stable platform.
- Extremely sturdy arm and legrests, wider footplates to create additional room and comfort.
- Corpus seating system for the ultimate in comfort with unmatched power tilt range of 0 - 45° and unprecedented recline of 85 - 150°.
- Ultra-strong actuators allow the user to safely elevate their feet above their hearts, while durable mesh-style upholstery helps keep moisture away from the body.



Specifications	
SWL	204kg
Seat Depth	19" - 25" (by 1" incr.)
Seat Width	19", 21", 23", 25"

## Evoke 2 HD

Code: WCFOLDHD 

- Designed for the larger user who wants a lightweight, highly adjustable manual wheelchair.
- Folding Frame.
- Highly adjustable, lightweight, aluminium wheelchair.
- Comprehensive range of compatible accessories is available for increased independence and enhanced clinical outcomes.



Specifications	
SWL	180kg
Seat Width	22" - 24"
Seat Depth	16" - 19"

# Wheelchair Range

## Minimaxx

Code: WCMINIMAX RENTAL HD

- Large rear wheels, with solid rubber tyres and hand rims, designed to give patients the option of manoeuvring and propelling themselves.
- Folding frame.
- Each rear wheel fitted with hand operated wheel brake to secure the Minimaxx when parked.
- Footrests designed to carry 100kg each.



Specifications	
SWL	325kg
Seat Width	61cm
Seat Depth	48cm

## Tilt in Space Junior

Code: WCTILT JNR 

The Rehab RX Junior is designed for the increased postural support and pressure care of children, via easy repositioning and advanced adjustability.

- Depth and width adjustable seat sizing, height adjustable, removable armrests, adjustable headrest.
- Tilt (26°) and backrest recline (30°).
- Contoured, built-up backrest and seat cushion.
- Advanced articulating leg rests closely mimic normal movement at the knee upon elevation.



Specifications	
SWL	180kg (Paediatric)
Seat Width	13" - 15"
Seat Depth	13" - 15"

## Wheelchair Accessories

The following are available for most wheelchairs. We have many options so please discuss these with us.

Code: - ELEVATING-LEG-WC RENTAL Elevating Leg Rest 

Code: - GUTTER-ARM-WC RENTAL Gutter Arms 

Code: - STUMP-SUPPORT-WC RENTAL Stump Support 

# Cushion Range

Permobil proudly supplies the ROHO® range of cushions, which are designed to move with the body and constantly adapt to the changes that happen throughout each day and over time.

## ROHO® Bariatric Cushion

Code: CUSHION ROHO BARI 

Sizes: 12 (22"), 13 (24"), 14 (26"), 15 (28")

- QUADTRO SELECT is the most clinically advanced cushion in the ROHO® product family, constantly adapting to needs or activities. The exclusive ISOFLO Memory Control and independent four-compartment design provides the ultimate control over the inflation levels across the cushion, making it easy to ensure proper positioning.
- The interconnected air cells envelop the body, delivering maximum immersion and positioning, while protecting from skin and soft tissue breakdown.



### LOW PROFILE

The interconnected air cells deliver sitting stability. It could be used as a secondary cushion.



### MID PROFILE

The mid-sized interconnected air cells deliver the ideal combination of customisable positioning and stability.



### HIGH PROFILE

The tall, interconnected air cells envelop your body, delivering maximum immersion for the ultimate in skin protection, positioning, stability and comfort.

While we have referred to a specific bariatric option above, we have a number of cushion and pressure support options that may also be suitable. Contact our rental team to discuss what is available and would suit your needs.

# ROHO® FAQs

## How do I know if a ROHO® is correctly inflated?

If you received a user manual with the ROHO® product, please refer to this. There may also be instructions inside the box your cushion or mattress came in.


If you did not receive a manual, do a search on our website for the relevant manual.  
[www.permobil.com/en-nz](http://www.permobil.com/en-nz)

## The ROHO® keeps going flat – what should we do?

Firstly, check the valve(s) on the ROHO® cushion – is this closed? If the valve is closed (and is remaining closed during the day) the ROHO® may have a puncture. We recommend that you seek assistance from the clinician or supplier who supplied you with the cushion as soon as possible.

## When I inflate the cushion, all the air remains in one corner of the cushion – what is going on?

This cushion is likely to be a ROHO® Quadro Select cushion. This cushion has an isoflo valve to limit air movement through the cushion. This cushion may need to be "reset". If you are unsure how to do this, please contact the cushion issuer.

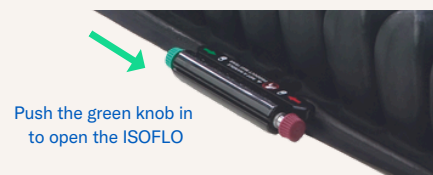

## The cushion and cover are soiled – can we wash them?

Yes! Please see the cleaning instructions or user manual provided. Alternatively, use this QR code to access our ROHO® Skills page on our website.



## What covers can we use on the ROHO® cushion or mattress section?

For cushions, please use the cover provided only. If keeping the cover clean is an issue, please request a second cover.

For mattress sections, please take care with bed linen used. Ideally, use a sheet with some stretch or a loosely fitted sheet to allow the person to immerse into the mattress section. Use of incontinence sheets is to be avoided if possible, as these will reduce the pressure relieving qualities of the cushion or mattress section.

## What is the weight limit on a ROHO® cushion or mattress overlay?

There is no weight limit on a ROHO® cushion or mattress, if appropriately sized.

# Chair Range

Our chair range provides additional comfort and support for a wide variety of users.

## MAXI Deluxe V2 Cloud Chair

Code: CLOUD BARI (ELECTRIC) 

The Maxi Deluxe V2 Cloud Chair offers ease of handling for carers and maximum comfort for patients.

- Power assisted tilt of backrest, seat and leg rests.
- Adjustable seat (tilt in space).
- Adjustable back, upright seated to full recline sleep position.
- Adjustable wings and drop-down arms for easy transfer.
- Adjustable leg support with optional foot rest.
- Locking castors - 3 x standard / 1 x directional.



Specifications	
SWL	300kg
Seat Size	W 72cm x D 55cm
Overall Size	W 95cm x D 95cm x H 138.5cm

## Bariatric Electric Rise and Recline

Code: LAZYBOY BARI

- Ideal solution for heavier weight users, with progressive conditions, who require pressure care management, who have Oedema legs or fluid build-up in legs, who have difficulty in stand to walk position.



Specifications	
SWL	255kg
Size	Dimensions (W x D x H cm)
Large	71cm x 56cm x 46cm
X-Large	76cm x 56cm x 46cm

# Chair Range

## Maxi Day Chair

Code: MAXI BARI 

- Seat with Anterior Tilt. Facilitates standing from sitting and provides greater comfort for stiff or painful hips.
- Height Adjustable. Easy push button adjustment for comfortable positioning and changing needs. Depth and backrest height adjustable with ease.
- Push handle and transport wheels. Allows a carer to tilt the chair backwards on the rear leg wheels for ease of transport (without a patient in the chair).



Specifications	
SWL	310kg
Seat Width	60cm - 81cm
Seat Depth	39cm - 53cm

## XXL Rehab Dining Chair

Code: BARI DINING RENTAL 

- The arm rests bend out like a butterfly, so body is not squeezed when getting in/out of the chair.
- Dining chair also provides better functionality for the user, compared to other chairs with straight and short armrests where it is almost impossible for a bariatric user to get a good grip.



Specifications	
SWL	325kg
Seat Size	W 61cm/71cm x D 40cm - 52cm
Seat Height	44cm-59cm

## Support Chair - Bariatric Chair Multiline XXL

Code: BARI TXCHAIR 

- Useful for patient transport, examination and treatment and post-op recovery.
- Capable of multiple positions, including lay flat position.
- All adjustments electronically controlled.



Specifications	
SWL	300kg
Seat Size	W 70cm x D 57cm
Overall Width	89cm

# Bathroom Range

We have a variety of options to support users in bathroom or wet areas.

## MAXI Shower Commode

Code: MAXI 65 SP Commode **P**

Durable and long lasting with a stainless steel frame, impervious PVC straps and water resistant castors.

- Simple to use swing up arms with padded arm rests and adjustable height swing away foot rests.
- Back rest with three flexible straps for comfort.
- Integrated push handle.
- 12.5cm braking castors.
- Durable and long lasting with a stainless steel frame, impervious PVC straps and water resistant castors.



Specifications	
SWL	300kg
Seat Size	W 65cm x D 48cm

## Bariatric Shower Commode

Code: BARI COMMODOE **P**

- Super heavy duty bariatric design.
- Anti-bacterial frame coating for improved infection control.
- Removable arms to assist with positioning and allow lateral slide transfers.
- Height adjustable, swing away footrests for optimal positioning and ease of transfers.
- Padded seat and arms for maximum comfort.



Specifications	
SWL	350kg
Seat Size	W 79cm x D 55cm
Overall Width	90cm

# Bathroom Range

## Weight Bearing Platform Shower Commode


Code: BARI SHOWERCOM 

- Heavy duty reinforced wide frame.
- 4 x total locking stainless steel castors provide brake access from all angles.
- Height adjustable reinforced armrests offer drop flat positioning to assist with lateral slide transfers.
- Weight bearing platform.
- The footrest slides completely away for transfers.



Specifications	
SWL	400kg
Seat Size	W 60cm x D 48cm
Overall Width	75cm

## Tilt and Recline Shower Commode

Code: SHOWERCOM TILT 

- 30° of backrest angle adjustment in 10° increments for effortless repositioning and maximum user comfort.
- -5° to 35° of gas-assisted tilt-in-space functionality minimises manual handling efforts during showering process.
- 4x total locking stainless steel castors provide brake access from all angles.



Specifications	
SWL	250kg
Seat Size	W 46cm x D 46cm
Overall Width	63cm

# Bathroom Range

## Shower Commode Maxi Tilt (Electric)

Code: SHOWERCOM ETILT 

- Assists transportation and positioning of the larger user for showering and can be positioned over the toilet for toileting.
- Electric tilt-in-space mechanism allows the carer to adjust the vertical position from 0° to 30° ensuring user comfortably retains correct seating position.
- Large calf pads to support the legs.
- Large backrest with headrest extension to support the upper body and head.



Specifications	
SWL	225kg
Seat Size	W 60cm x D 47cm
Overall Width	71cm

## XXL Rehab Shower Bench

Code: BARI SHOWERBENCH 

- Strongly engineered and functionally designed exclusively for bariatric users.
- The height of the chair is adjustable from 44cm to 59cm.
- Also available as a commode with pan.



Specifications	
SWL	325kg
Seat Size	W 61cm x D 40-52cm
Overall Width	69cm

# Manual Handling Range

Our manual handling range supports confident and comfortable movement and transfers.

## MAXI Walking Tutor Frame

Code: MAXI WALKER 

- SWL - 250kg.
- Rigid, fully welded, four wheeled, mobile walking frame with adjustable height forearm supports with upright foam padded handles.
- Supports also adjust in cuff length and angle.
- Frame has swivel front castors and rear wheels or optional glides.
- Padded chest and arm rest is available as an option.



## Hoist (200kg) - Standing Patient Lifter


Code: HOISTA200S RENTAL 

- SWL - 200kg.
- Compact, lightweight and adjustable standing lifter ideal for home, community and facility environments.



# Manual Handling Range

## Hoist (320kg) - Aluminium Bariatric Lifter

Code: HOISTA320 RENTAL 

- SWL - 320kg.
- The broad base and ample patient clearance compliments the substantial safe working load.



## Slings for Hoists

Code: SLING 

- Our range includes general purpose and bariatric slings with or without head support, hygiene slings for toileting and dressing, all-day slings for those who need to remain in their sling for extended periods and standing slings for patients who can bear weight.
- Size: XS, S, M, L, XL, XXL




# Permobil Clinical Resources

Our rental team work closely with clients and mobility experts to ensure that suggested products will meet the identified needs and genuinely support the user effectively. We recommend that you speak to your clinician or therapist to further understand the benefits of products and best use of them.

The below clinical resources are free to access for all healthcare professionals and consumers. Simply scan the QR codes to download.

This purpose of these guides are to:

- Foster a more holistic approach to the management and prevention of pressure injuries (PI) by simultaneously treating supine and seated postures.
- Review and adopt the International Clinical Practice Guideline (CPG) as published by the Pan Pacific Pressure Injury Alliance (PPPIA).
- Adopt key recommendations from the PPPIA's CPG when considering support surfaces and/or seating components.
- Understand the critical role wheelchair configuration and equipment choices play in the prevention and treatment of pressure injuries.
- Provide brief overviews, best used in conjunction with the full 2019 International Clinical Practice Guideline and the free abridged Quick Reference Guide.

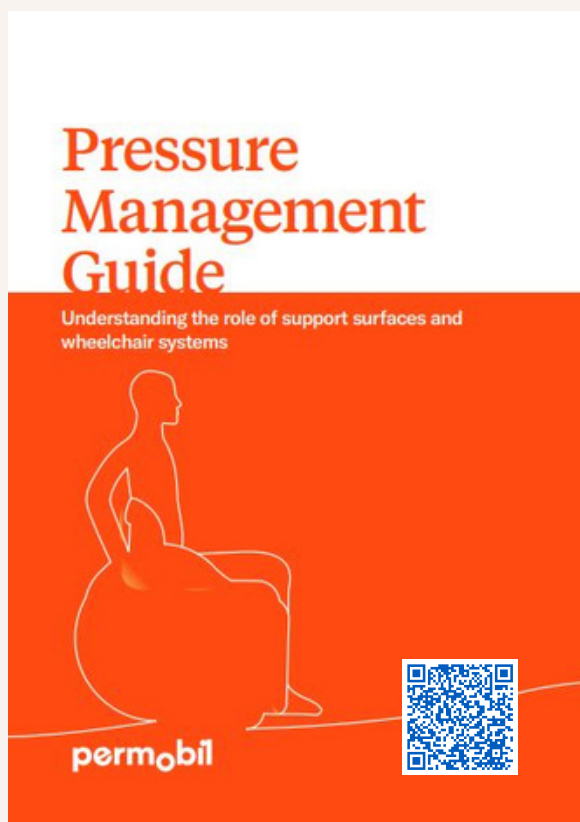

# FAQs and Contacts

## Which areas of New Zealand do you support?

We provide a nationwide service. Our three branches are based in Auckland, Wellington and Christchurch. We deliver within a large radius of those sites and can ship to all regions. Call our team to discuss options.

## What if I get the item and it's not quite right?

This happens from time to time – simply call our team, explain the issue and they'll work with you to identify a solution.

## What are your charges?

Once our team confirm the product(s), charges will be set out and discussed. These include a daily rental rate, alongside delivery and collection charges (where applicable). A quote is provided before confirmation of the hire. There are various payment options depending on needs, which the team can discuss with you.

## What if we need training or further product information?

Discuss this with the team when you are in contact – we provide training to ensure users receive maximum benefit from rental products.

**Pharmac:** Items marked with **P** are listed on Pharmac contract PML-2019-12-12

## Contact Us

### Auckland Rental

Rental.NZ@permobil.com

0800 115 222

181A McLeod Road, Te Atatu South

Auckland 0610, New Zealand

### Wellington Rental

RentalNZ.Nth@permobil.com

0800 363 123

15 Raroa Road, Hutt Central

Lower Hutt 5010, New Zealand

### Christchurch Rental

RentalNZ.Sth@permobil.com

03 354 9239

3/114 Sawyers Arms Road, Northcote

Christchurch 8052, New Zealand

### Permobil Contract Manager

Susie Ballantyne

Commercial Services Director

Susie.Ballantyne@permobil.com

021 246 1885



permobil

[www.permobil.com/en-nz](http://www.permobil.com/en-nz)

BR-NZ-BARI-CAT-01.0